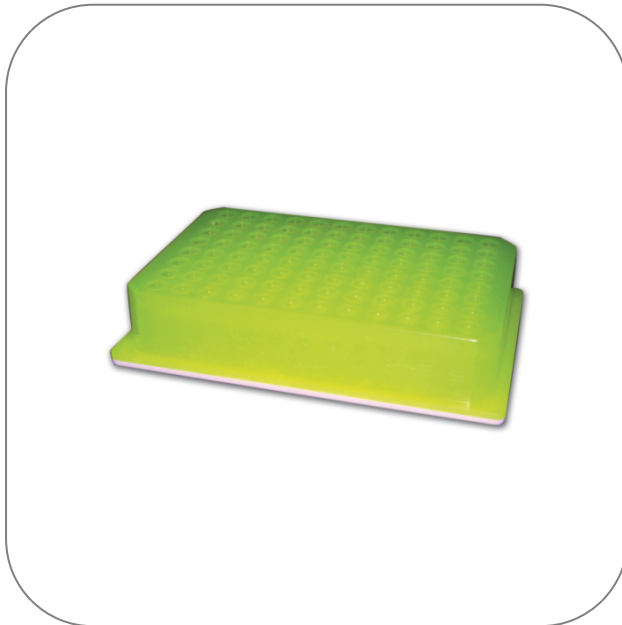




## StarChill SBS, Pink



### Product Highlights

- 96-well rack (no lid), for use on AUTOMATED SYSTEMS with: 0.2 ml tubes (single or strip) or 0.2 ml 96-well PCR Plate
- Temperature colour change: 7 °C
- Estimated time at 4 °C: 2+ hours\*
- Dimensions (WxDxH): 128 x 85 x 26 mm
- StarChill Rack must be frozen for a minimum of 8 hours at -20°C before use
- Freeze the StarChill Rack upside down in order for the gel to freeze evenly

## Product Applications

Sample storage

STARLAB reserves the right to make changes at any time and without prior notice. The content and design of this PDF are protected by national and international copyright law and are the property of STARLAB International GmbH. Any duplication, editing, distribution and any kind of use and utilization of this PDF content in electronic systems, online media and / or libraries or similar databases requires the prior consent of STARLAB International GmbH.

**Starlab (UK) LTD**  
5 Tanners Drive  
MK14 5BU Milton Keynes  
Email: [info@starlab.co.uk](mailto:info@starlab.co.uk)



## General Data

Art. No.	E2396-5654
Pack Size	2 Pcs. (1 Pack × 2 Pcs.)
Dimensions (W × D × H)	128 mm × 85 mm × 26 mm
Dimensions (W × D)	128 mm × 85 mm
Autoclavable	No
Number of wells	96 wells
Rack Type	Chiller Rack
Color	Pink (warm)/purple (cold)

STARLAB reserves the right to make changes at any time and without prior notice. The content and design of this PDF are protected by national and international copyright law and are the property of STARLAB International GmbH. Any duplication, editing, distribution and any kind of use and utilization of this PDF content in electronic systems, online media and / or libraries or similar databases requires the prior consent of STARLAB International GmbH.

**Starlab (UK) LTD**  
5 Tanners Drive  
MK14 5BU Milton Keynes  
Email: [info@starlab.co.uk](mailto:info@starlab.co.uk)



## More informations about StarChill SBS, Pink

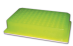

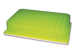

### StarChill high performance chiller racks

Maintain a fridge environment during automation! The centre section of StarChill racks changes colour when frozen. It then gradually turns back to its original colour as it warms up, thus giving an instant, visual-indication of the temperature.

\*Time will vary, depending on ambient temperature. All StarChill racks must be frozen for a minimum of 8 hours at -20 °C before use. Not Autoclavable.



## All Variations

PRODUCT NAME		PACKAGING SIZE	ART. NO.
	<a href="#">StarChill SBS, Pink</a> <b>Color:</b> Pink (warm)/purple (cold)	2 Pcs. (1 Pack × 2 Pcs.)	E2396-5654
	 <a href="#">Create PDF Data Sheet</a>		
	<a href="#">StarChill SBS, Green</a> <b>Color:</b> Yellow (warm)/green (cold)	2 Pcs. (1 Pack × 2 Pcs.)	E2396-5656
	 <a href="#">Create PDF Data Sheet</a>		

STARLAB reserves the right to make changes at any time and without prior notice. The content and design of this PDF are protected by national and international copyright law and are the property of STARLAB International GmbH. Any duplication, editing, distribution and any kind of use and utilization of this PDF content in electronic systems, online media and / or libraries or similar databases requires the prior consent of STARLAB International GmbH.

**Starlab (UK) LTD**  
5 Tanners Drive  
MK14 5BU Milton Keynes  
Email: [info@starlab.co.uk](mailto:info@starlab.co.uk)