

New!

STARLAB *focus*

TipOne® Repeat Dispenser Tips

You want more freedom? No problem! Our new TipOne® Repeat Dispenser Tips are available in an extra volume of 0.1 ml! Additionally our new TipOne® Repeat Dispenser Tips fit more dispenser models than before (e. g. even now with the Multipette® Plus). This new range is still made from the same high-quality materials, but now have the following features:



What's new?

- /// Compatible with a wider range of dispensers
- /// High accuracy and precision
- /// Tip volume recognition with Multipette® Plus
- /// Now available in eight volumes – new 0.1 ml tip!
- /// Up to 100 dispenser steps with Multipette® Plus
- /// Separate adapters for 25 ml and 50 ml tips, easily distinguishable by colour
- /// Extended piston length of most tips minimises dead air volume, improves reagent delivery and provides easier access to narrow vessels
- /// Lower dispensing forces

TipOne Repeat® Compatibility

STARLAB	RepeatOne® Dispenser
Brand	Handystep® Manual, HandyStep® S, HandyStep® electronic
Denville	Faststepper®
Eppendorf	Multipette® 4780 und Multipette® Plus*, Multipette® stream/Xstream**
Gilson	Pipetman®, Distriman®
Rainin	AutoRep® E, AutoRep® S, AutoRep® M
Ritter	Ripette®, Ripette® pro, Ripette® genX
Sartorius	Biohit® Mechanical Stepper, Biohit® eLine® Pro, Biohit® eLine® Lite
Scilogex	Stepmate®
Thermo Fisher	Fisherbrand® Repeat Stepper, FinnpiPETTE® Stepper

* Tip recognition works only with old sensor (prior to service).

**For use with models manufactured prior to Dec 2013 only.

STARLAB International GmbH

Neuer Höltigbaum 38
22143 Hamburg
Germany

T: +49 (0)40 675 99 39 0
F: +49 (0)40 675 99 39 20
info@starlab.de

Passionate for Science.



www.starlabgroup.com

New!


Order Information

Dispenser Tips		
Cat. No.	Description	Pack Size
Non-Sterile		
S4751-0010	0.1 ml TipOne® Repeat Dispenser Tips	100
S4751-0050	0.5 ml TipOne® Repeat Dispenser Tips	100
S4751-0100	1.0 ml TipOne® Repeat Dispenser Tips	100
S4751-0250	2.5 ml TipOne® Repeat Dispenser Tips	100
S4751-0500	5.0 ml TipOne® Repeat Dispenser Tips	100
S4751-1000	10 ml TipOne® Repeat Dispenser Tips	100
S4751-2500	25 ml TipOne® Repeat Dispenser Tips*	100
S4751-5000	50 ml TipOne® Repeat Dispenser Tips*	100
Sterile, Individually Wrapped		
S4761-0010	0.1 ml TipOne® Repeat Dispenser Tips (Sterile)	100
S4761-0050	0.5 ml TipOne® Repeat Dispenser Tips (Sterile)	100
S4761-0100	1.0 ml TipOne® Repeat Dispenser Tips (Sterile)	100
S4761-0250	2.5 ml TipOne® Repeat Dispenser Tips (Sterile)	100
S4761-0500	5.0 ml TipOne® Repeat Dispenser Tips (Sterile)	100
S4761-1000	10 ml TipOne® Repeat Dispenser Tips (Sterile)	100
S4761-2500	25 ml TipOne® Repeat Dispenser Tips (Sterile)*	100
S4761-5000	50 ml TipOne® Repeat Dispenser Tips (Sterile)*	100
Adapters		
S4756-2500	Adapter 25 ml	1
S4756-5000	Adapter 50 ml	1
S4765-2500	Adapter 25 ml (Sterile)	7
S4765-5000	Adapter 50 ml (Sterile)	7

*25ml and 50ml tips are supplied with four adapters.
NOTE: TipOne® Repeat dispenser tips are not autoclavable.

Repeat Dispenser

RepeatOne® allows the selection of five separate volumes from each TipOne® Repeat Dispenser Tip.

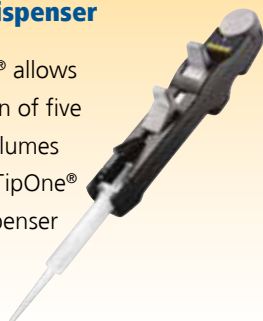

Features

- /// Available non-sterile bulk, or sterile individually wrapped
- /// High accuracy and precision
- /// ISO 8655 manufactured
- /// Polypropylene cylinder
- /// Polyethylene piston
- /// All TipOne® Repeat tips are certified RNase, DNase, DNA & Pyrogen free

Application Guide

Dispenser Tip Volume	Dispenser Adjustment				
	1	2	3	4	5
0.1 ml	2 µl	4 µl	6 µl	8 µl	10 µl
0.5 ml	10 µl	20 µl	30 µl	40 µl	50 µl
1.0 ml	20 µl	40 µl	60 µl	80 µl	100 µl
2.5 ml	50 µl	100 µl	150 µl	200 µl	250 µl
5.0 ml	100 µl	200 µl	300 µl	400 µl	500 µl
10.0 ml	200 µl	400 µl	600 µl	800 µl	1000 µl
25.0 ml	500 µl	1000 µl	1500 µl	2000 µl	2500 µl
50.0 ml	1000 µl	2000 µl	3000 µl	4000 µl	5000 µl
Steps*	48	23	15	11	8

*Number of steps may vary according to dispenser manufacturer.

STARLAB International GmbH
Neuer Höltigbaum 38
22143 Hamburg
Germany

T: +49 (0)40 675 99 39 0
F: +49 (0)40 675 99 39 20
info@starlab.de

Passionate for Science.



www.starlabgroup.com